

EASY PULSE


Efficient and compact mechanical CPR



EASY PULSE

Improving CPR quality with SCHILLER's EASY PULSE:

- ∴ Ideal for confined spaces and aircraft, thanks to its small size and light weight
- ∴ May be used in a head-up position of up to 45°
- ∴ Delivers more consistent rate and depth than manual chest compressions

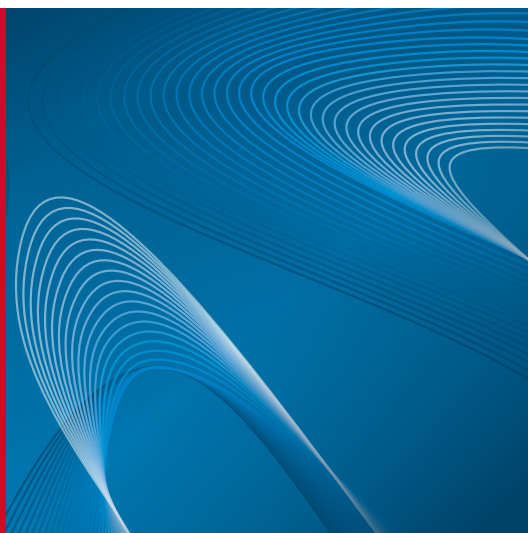

EASY PULSE provides efficient and consistent CPR during the rescue process

IMPROVE CPR QUALITY AND EFFICIENCY

Performing manual chest compressions well for an extended period of time is very challenging. With EASY PULSE, emergency responders can focus on other critical tasks while the device provides automatic chest compressions at a consistent rate and depth. SCHILLER's EASY PULSE is a portable, stand-alone solution for more effective resuscitation.

CPR IN ANY SITUATION

EASY PULSE is directly attached to the patient's upper body and can thus be used in any situation. Thanks to its small size, EASY PULSE is also ideal for confined spaces in disaster situations and accident sites. CPR can be performed even during patient rescue.

CPR IN HEAD-UP POSITION

EASY PULSE can be used in a head-up position of up to 45°, ensuring consistent CPR quality even during stair or rope rescues.

IDEAL FOR AIRCRAFT

EASY PULSE is such a simple device that a single person can operate it. Its size and weight are unparalleled: EASY PULSE opens up entirely new possibilities, e.g. for air rescue services.

SOPHISTICATED TECHNOLOGY

The EASY PULSE has a unique combination of piston and band which allows a multidirectional 3D compression and therefore maximum efficiency.

The weight of the device is evenly distributed over the baseplate and not on a single point and therefore has no influence on chest compressions or on the recoil.

TWO TYPES OF SLIDER



The double slider complements the standard version by providing a wider base for enhanced stability during extended periods of use. Two additional safety straps complement the proven fastening system and ensure secure positioning even under the most demanding operating conditions.

CHARGER

One battery can be charged in the left slot in fast charging mode (80% charged in 45 minutes). Two batteries can be charged simultaneously.





POWER SUPPLY

Replaceable battery with charge level indicator. Autonomous operation for 45 minutes (with one battery). Option to change the battery during operation when power is supplied by the vehicle or mains.



EASE OF USE

Pressing two buttons is enough to initiate resuscitation immediately

SMALL SIZE

Ideal for limited spaces, such as in disaster relief operations

POWER CONNECTION CABLE

Power supply for 100VAC to 240VAC input voltage. EASY PULSE can be directly operated using a 12VDC to 48VDC power source via a connection cable.

LIGHTWEIGHT

Unbelievably small and lightweight CPR system



FAST HOOK-UP


Colour-coded straps to speed up patient hook-up. It only takes 15 seconds for trained personnel to attach the device thanks to the convenient buckle system.

EASY TO MOUNT

Thanks to a slider and buckle system, it can be easily attached to the upper body. Easy to position and intuitive to use.

COMPREHENSIVE

Option of 30:2 compression-ventilation cycle

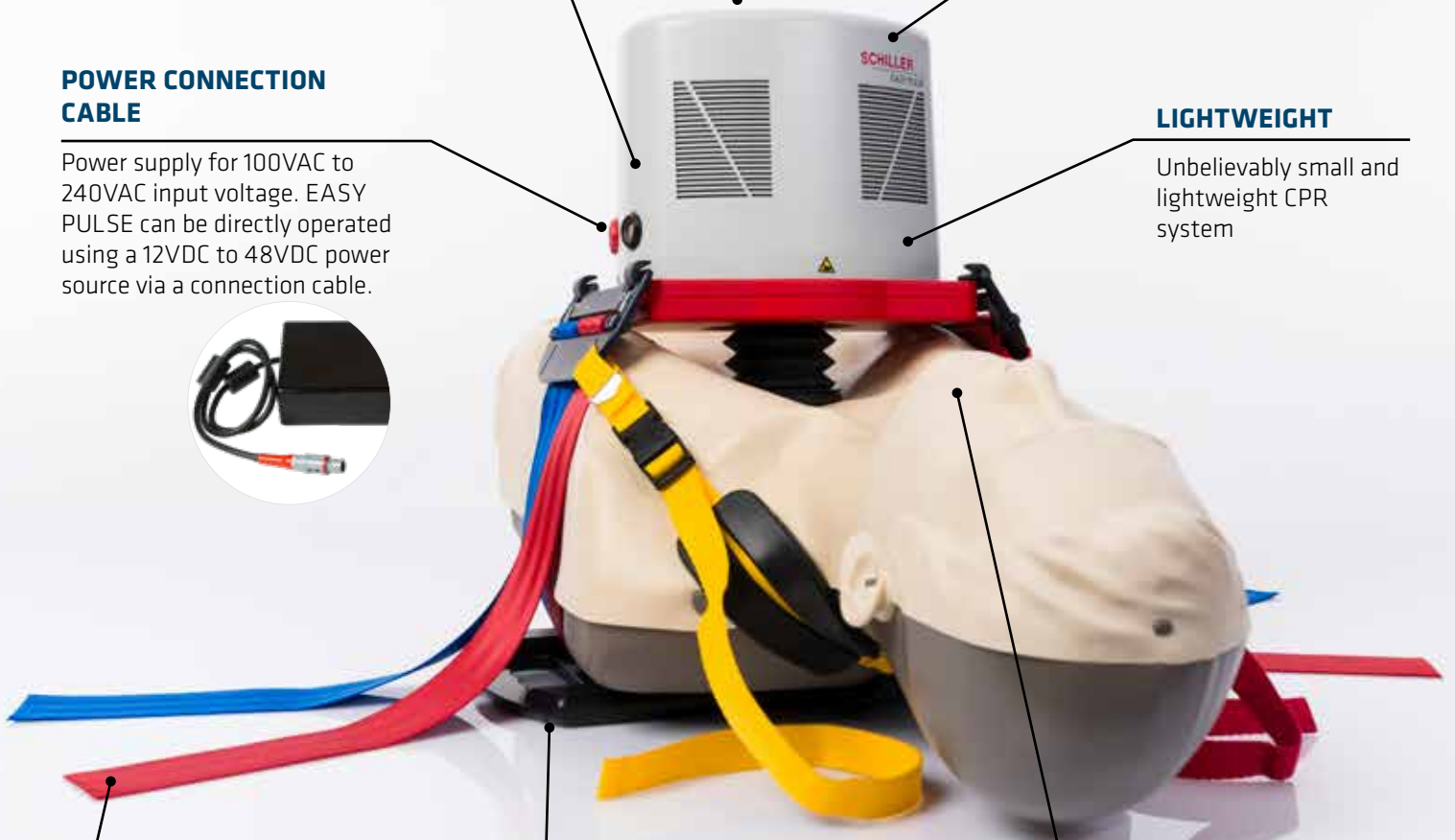



Americas
SCHILLER Americas Inc.
 Doral, Florida 33172

North America:
 Phone +1 786 845 06 20
 Fax +1 786 845 06 02
sales@schilleramericas.com
www.schilleramericas.com

Latin America & Caribbean
 Phone +1 305 591 11 21
 Fax +1 786 845 06 02
sales@schilleramericas.com
www.schilleramericas.com



Asia
SCHILLER Asia-Pacific / Malaysia
 52200 Kuala Lumpur, Malaysia
 Phone +603 6272 3033
sales@schiller.com.my
www.schiller-asia.com



Austria
SCHILLER Handelsgesellschaft m.b.H.
 A-4040 Linz
 Phone +43 732 70 99 0
sales@schiller.at
www.schiller.at



China
Alfred Schiller (Beijing) Medical Technology Co., Ltd.
 100102 Beijing, China
 Phone +86 10 84565453
info@schillerchina.com
www.schillermedical.cn



Croatia
SCHILLER d.o.o.
 10000 Zagreb
 Phone +385 1 4059 559
info@schiller.hr
www.schiller.hr



France
SCHILLER Médical
 F-67160 Wissembourg
 Phone +33 3 88 63 36 00
info@schiller.fr
www.schiller-medical.fr



France (distribution France)
SCHILLER France S.A.S.
 F-77608 Bussy St Georges
 Phone +33 1 64 66 50 00
contact@schillerfrance.fr
www.schiller-france.com



Germany
SCHILLER Medizintechnik GmbH
 D-85622 Feldkirchen b. München
 Phone +49 89 62 99 81 0
info@schillermed.de
www.schillermed.de



India
SCHILLER Healthcare India Pvt. Ltd.
 Mumbai – 400 059, India
 Phone +91 22 6152 3333 / 2920 9141
sales@schillerindia.com
www.schillerindia.com



Poland
SCHILLER Poland Sp. z o.o.
 PL-02-729 Warszawa
 Phone +48 22 843 20 89 / +48 22 647 35 90
schiller@schiller.pl
www.schiller.pl



Russia & C.I.S.
AO SCHILLER.RU
 119049 Moscow, Russia
 Phone +7 (495) 970 11 33
mail@schiller.ru
www.schiller.ru



Serbia
SCHILLER d.o.o.
 11010 Beograd
 Phone +381 11 39 79 508
info@schiller.rs
www.schiller.rs



Slovenia
SCHILLER d.o.o.
 2310 Slovenska Bistrica
 Phone +386 2 843 00 56
info@schiller.si
www.schiller.si



Spain
SCHILLER ESPAÑA, S.A.
 ES-28232 Las Rozas/Madrid
 Phone +34 91 713 01 76
schiller@schiller.es
www.schiller.es



Switzerland
SCHILLER Schweiz AG
 CH-8912 Obfelden
 Phone +41 44 744 30 00
sales@schiller-schweiz.ch
www.schiller-schweiz.ch



Turkey
SCHILLER TÜRKİYE
 Okmeydanı-Sisli – Istanbul
 Phone +90 212 210 8681
info@schiller.com.tr
www.schiller-turkiye.com



United Kingdom
SCHILLER UK Ltd.
 Bromborough, Wirral,
 Merseyside, CH62 3RJ
 Phone +44 1333 312150
sales@schilleruk.com
service@schilleruk.com
www.schilleruk.com

SCHILLER reserves the right to make changes in design, specifications and features shown or discontinue the product described at any time without notice or obligation. Images may be representative. Some products may not be regulatory cleared or released for sale in all markets. Contact your SCHILLER representative for the most current information and availability.

